

**MICROGEN BIOPRODUCTS LTD**

---



## **PATH-CHEK™ HYGIENE PROTEIN SWAB**

**SIMPLE AND RAPID TEST FOR THE  
DETERMINATION OF THE LEVEL OF  
SURFACE PROTEIN IN FOOD HANDLING  
ENVIRONMENTS**

**Protecting Food and Health**

---

MICROGEN BIOPRODUCTS LTD

1 Admiralty Way  
Camberley  
Surrey  
GU15 3DT  
United Kingdom

Phone: +44 (0)1276 600081  
Fax: +44 (0)1276 600151  
E-mail: [sales@microgenbioproducts.com](mailto:sales@microgenbioproducts.com)  
[www.microgenbioproducts.com](http://www.microgenbioproducts.com)




The increased awareness of the need for high levels of hygiene in food handling and processing environments combined with the widespread implementation of HACCP Programs has resulted in the need for simple and effective methods for the rapid assessment of hygiene and the effectiveness of cleaning programs.

Pathchek™ Hygiene is a chemical test based on the detection of protein, a common component of most food types. This test is not targeted specifically at bacteria. When cleaning procedures are not effective or have been performed poorly, both food residues and bacteria will be present and will result in a positive reaction.

Pathchek™ Hygiene is provided in a simple peel open plastic pouch which contains a pre-moistened swab for sampling test surfaces. After sampling, the swab is returned to the pouch and pushed down into the developing reagent. If positive (indicating the presence of food residues and bacteria), a colour reaction (green to blue) will develop within 5 seconds.

### Path Chek Procedure:

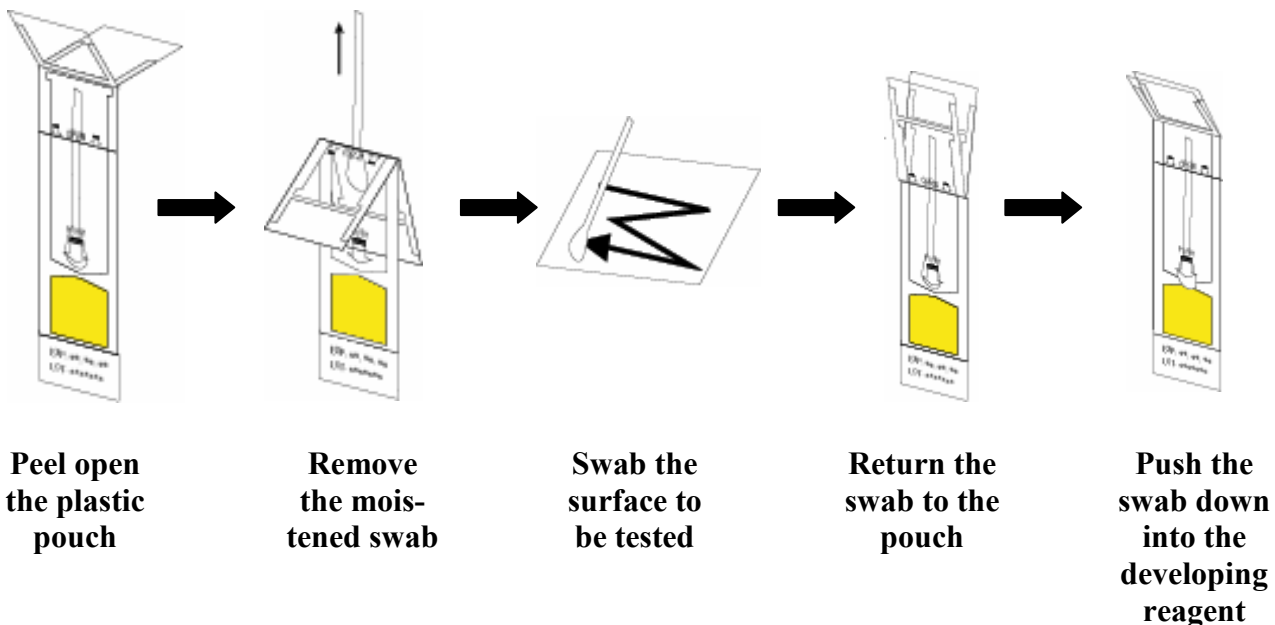

<b>Code</b>	PC005-100
<b>Description</b>	Pathchek™ Hygiene
<b>Size</b>	100 tests
<b>Storage</b>	Below 30°C